

Mobile IV Stand — instructions for use

Model: IVSTND-260 | English | Rev 1.2 | 25 October 2026

Read these instructions before assembly and use. Keep them with the stand. This is a non-sterile, reusable medical device for use by trained healthcare staff indoors.

Intended purpose

IVSTND-260 supports compatible infusion bags or bottles on four hooks and organises light clinical supplies on a removable tray and lower basket. It also has supports for a waste bin and a compatible 1 L round sharps container. Its contribution is mechanical support and access to supplies during care.

The stand does not pump, measure, administer or monitor an infusion and does not contact the fluid pathway. It does not support or transport the patient.

Operating limits

Feature	Limit or specification
Maximum combined load	10 kg across the hooks, tray and basket
All four hooks combined	4 kg total, distributed evenly
Removable tray	2 kg
Lower basket	4 kg
Height adjustment	Minimum 1,200 mm · normal 1,750 mm · maximum 2,000 mm Measured from the floor to the top of the hooks
Base	Five wheels, two fitted with brakes
Sharps container	Compatible round 1 L container; follow its own instructions

The hook, tray and basket ratings may be used together within the 10 kg combined limit. Do not use internal test loads as operating ratings.

Do not use the stand for


Sitting, standing, lifting, walking support or patient transport; mounting pumps or other equipment on the tube; outdoor use; stairs; uneven or sloping floors. MRI suitability has not been established.

Stop use if the stand is unstable or a joint, height lock, wheel, brake or attachment is defective. Residual risks include tipping or falling loads, unintended movement, pinching at clamps, snagged tubing, sharps injury and cross-contamination from inadequate cleaning and disinfection. Keep within all load limits, use the brakes, secure loads, keep hands and lines clear, and follow the sharps-container and cleaning instructions.

Assembly


Follow the sequence below using the matching supplied parts and hardware. Do not substitute fasteners.

1 Prepare the telescopic pole

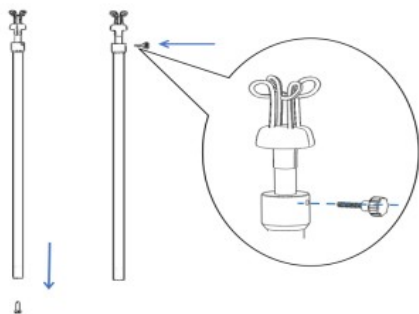


Remove the hexagonal screw from the bottom of the pole. Align the hand-tightening screw with the hole at the telescopic connection and tighten it to secure the pole. No tools required.

2 Fit the tray support

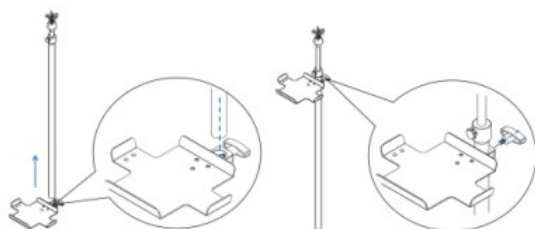


Slide the tray support onto the pole from the bottom and position it at the required height. Secure it with the black T-shaped screw supplied for the support. No tools required.

3 Attach the pole to the base

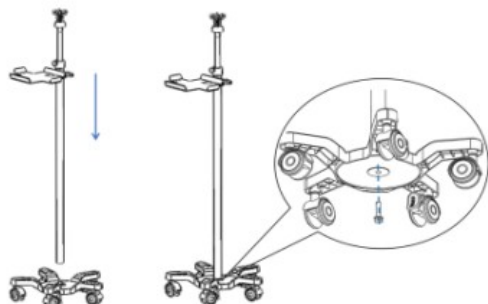

Insert the plastic cover into the central hole of the base. Insert the pole and secure it from underneath with the M10 × 35 hexagonal combination screw. Use the 16 mm open-ended wrench.


4 Fit the lower mounting plate

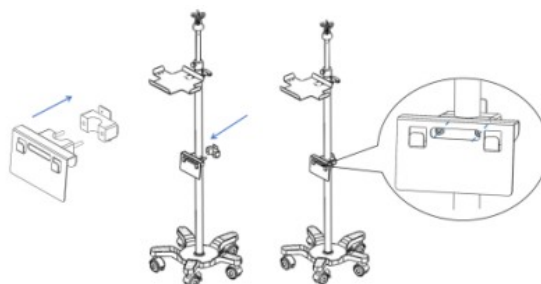

Remove the clamp from the lower mounting-plate assembly. Position the plate and clamp around the pole as shown and secure them with two M6 × 65 socket-head screws. Use the 5 mm hex key.


Assembly continued

5 Fit the basket

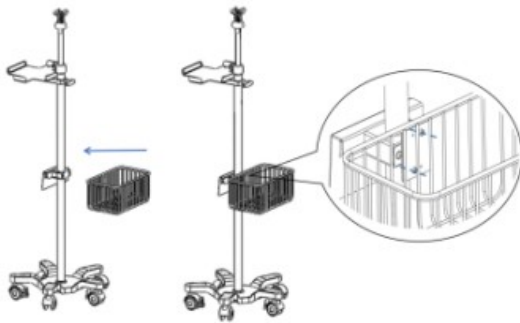


Align the basket with the threaded holes on the back of the mounting plate. Secure it with two M4 × 12 socket-head screws. Use the 3 mm hex key.

6 Fit the waste bin

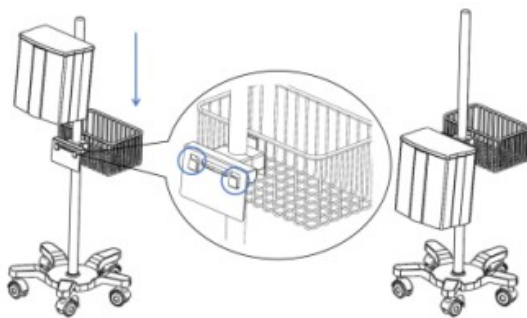


Align the waste bin with the two hooks on the mounting plate. Lower it onto the hooks and press down firmly. No tools required.

7 Fit the sharps-container holder

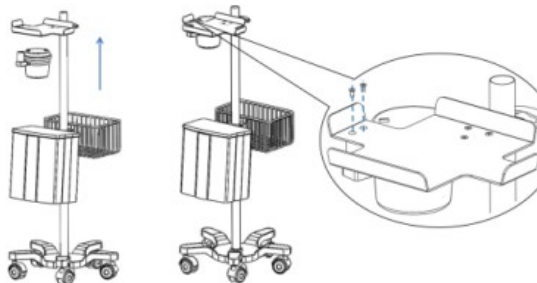

Position the holder beneath the tray support, aligning the two fixing holes shown. Secure it with two M4 × 10 countersunk hex-socket screws. Use the 2.5 mm hex key.

8 Fit the removable tray

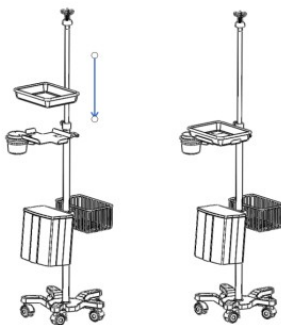

Lower the stainless-steel tray into the tray support. Check that it is seated securely before loading. No tools required.

After assembly, check every joint and attachment, both brakes and the height lock before placing the stand in service.

Before each use

- Check that the stand is clean, the identification is readable and the surfaces are intact.
- Check all joints, clamps and accessories for secure attachment. Check that the tray is seated properly.
- Check that all wheels turn freely and both brakes operate. Apply both brakes and check that the stand stays stationary.
- Support the upper tube before loosening the height lock. Adjust within 1200–2000 mm, retighten the lock and check that the tube does not slip before loading.
- Distribute the hook loads evenly and keep within all individual and combined load limits.

Stationary use and movement

Apply both wheel brakes during stationary use. Keep fingers clear of sliding and locking parts. Keep the infusion tubing clear of the wheels and any parts that could trap or pull it.

Before moving the stand, secure the containers and retract the sharps-container holder beneath the tray. Release the brakes and push slowly, keeping the load stable. Approach small door thresholds cautiously; stop and obtain assistance if the stand is unstable. Do not force the loaded stand over an obstruction.

Using the sharps-container holder

Pull out the holder for use and return it beneath the tray for movement. Use the compatible round 1 L sharps container. Follow the container manufacturer's fill limit, closure and disposal instructions. Do not place sharps in the ordinary waste bin. The holder does not replace the container's protective features or instructions.

Cleaning and disinfection

Remove supplies and waste safely. Remove soiling with a soft cloth and mild detergent in water. Wipe off detergent with a clean damp cloth and allow the surface to air dry. Do not mix chemicals, immerse the stand, pressure-wash it or use abrasive cleaners.

Disinfect with ready-to-use hospital surface-disinfectant wipes based on 70% isopropyl alcohol. Keep the surface visibly wet for the full contact time stated on the wipe label, re-wetting if needed. Do not use other disinfectant types.

Maintenance, service life and warranty

Specific scheduled maintenance: none. Carry out the before-use checks in this manual. Stop using a stand with loose joints, a slipping lock, faulty brakes, damaged parts or surfaces that cannot be cleaned. Qualified personnel must repair it and verify safe function before reuse.

Expected service life: at least five years under normal use, subject to the condition checks in this manual. Warranty: two years. The expected service life is not an extension of the warranty.

End of use

Decontaminate the stand before disposal. Handle sharps and clinical waste separately under the applicable local rules and the container instructions. Dispose of the stand's materials in accordance with local requirements.

Identification and reporting

Model: IVSTND-260. GTIN: 6980438430027. Identify an individual unit by its serial number and lot number on the product label. Quote the model, serial number and lot when contacting EverCare about a defect or incident. The 14-digit UDI-DI representation is 06980438430027. The house symbol on the label means "Indoor use only".

Report any serious incident associated with the device to the manufacturer and to the competent authority of the Member State where the user or patient is established.

Manufacturer

EVERCARE INNOVATIONS (SHENZHEN) CO., LTD.

Unit 102, Building J, No. 3 Xiang Yuan Pai Road, Pingxi Community, Pingdi Subdistrict, Longgang District, Shenzhen, Guangdong 518000, China

Telephone: +86 136 3167 8692

Email: service@evercare-innovations.com

Website: www.evercare-innovations.com

European authorised representative

MDSS GmbH

Schiffgraben 41, 30175 Hannover, Germany

Document ECI-IFU-STND-260 | Rev 1.2 | 25 October 2026

Identification and document control

Mobile IV Stand

Medical device | Non-sterile | Reusable | Professional healthcare use

Model / catalogue number	UDI-DI (GTIN, 14-digit representation)
IVSTND-260	06980438430027

The UDI-DI identifies the product model, not an individual unit. Read the unit's serial number and lot number from its product label. Match the model and unit identification on the device, packaging and accompanying records before use. Do not substitute the manufacturer SRN for a product UDI.

Manufacturer identification

EVERCARE INNOVATIONS (SHENZHEN) CO., LTD.

Unit 102, Building J, No. 3 Xiang Yuan Pai Road, Pingxi Community, Pingdi Subdistrict, Longgang District, Shenzhen, Guangdong 518000, China

Manufacturer SRN: CN-MF-000060301

Email: service@evercare-innovations.com | Telephone: +86 136 3167 8692

Website: www.evercare-innovations.com

EU authorised representative

MDSS GmbH

Schiffgraben 41, 30175 Hannover, Germany

Revision identification

ECI-IFU-STND-260 | Rev 1.2 | English | 25 October 2026

Changes since released Rev 1.0: full manufacturer address and identification page (manufacturer SRN, model and GTIN, authorised representative and CE marking); residual risks; and cleaning and disinfection instructions. The indoor-use-only label symbol is explained. Revision 1.1 was an unreleased draft.

